## RH 600

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## Contents

### Page

|  |           |
|--|-----------|
| <b>For Your Safety and That of Your Patients .....</b> | <b>31</b> |
| <b>Intended Use.....</b>                               | <b>32</b> |
| <b>Preparation.....</b>                                | <b>34</b> |
| Positioning the apparatus.....                         | 34        |
| Checking readiness for operation .....                 | 34        |
| <b>Operation .....</b>                                 | <b>37</b> |
| Safety precautions.....                                | 37        |
| Setting the heater output manually .....               | 38        |
| Using skin temperature control .....                   | 39        |
| Using the timer .....                                  | 43        |
| Adjusting the halogen lamp beam .....                  | 44        |
| Shutdown.....  | 44        |
| Setting instructions.....                              | 45        |
| Phototherapy .....                                     | 47        |
| <b>Care.....</b>                                       | <b>48</b> |
| <b>Maintenance intervals .....</b>                     | <b>49</b> |
| Disposal of used batteries.....                        | 49        |
| <b>Troubleshooting .....</b>                           | <b>50</b> |
| <b>What's what.....</b>                                | <b>51</b> |
| <b>Technical Data .....</b>                            | <b>53</b> |
| <b>Construction and description.....</b>               | <b>56</b> |
| <b>Order List .....</b>                                | <b>57</b> |
| <b>Parts List .....</b>                                | <b>57</b> |

## For Your Safety and That of Your Patients

### Strictly follow the Instructions for Use

Any use of the apparatus requires full understanding and strict observation of these instructions. The apparatus is only to be used for purposes specified here.

### Maintenance

The apparatus must be inspected and serviced yearly by trained service personnel (and a record kept). Repair and general overhaul of the apparatus may only be carried out by trained service personnel. We recommend that a service contract be obtained with DrägerService and that all repairs also be carried out by them. Only authentic Dräger spare parts may be used for maintenance. Observe chapter "Maintenance Intervals".

### Accessories

Do not use accessory parts other than those in the order list.

### Not for use in areas of explosion hazard

This apparatus is neither approved nor certified for use in areas where combustible or explosive gas mixtures are likely to occur.

### Liability for proper function or damage

The liability for the proper function of the apparatus is irrevocably transferred to the owner or operator to the extent that the apparatus is serviced or repaired by personnel not employed or authorized by DrägerService or if the apparatus is used in a manner not conforming to its intended use.

Dräger cannot be held responsible for damage caused by non-compliance with the recommendations given above. The warranty and liability provisions of the terms of sale and delivery of Dräger are likewise not modified by the recommendations given above.

Dräger Medical AG & Co. KGaA

## Intended Use

The RH 600 infant radiant heater is a radiant heater for keeping premature babies, neonates and infants, up to a body weight of 8 kg, warm.

It is intended for use in the delivery room, neonate unit, paediatric unit and paediatric intensive care units.

This apparatus must be operated by a doctor or by nurses or care assistants under the instructions of a doctor.

### Therapeutic and care possibilities:

- Reanimation station
- Care station for normal and intensive care
- Warming and compensating for heat loss
- Weaning of incubator infants
- Timer for APGAR test

### Heating possibilities of the RH 600 infant radiant heater:

- Operation of the radiant heater with set heat output:  
Radiant heater output controlled in manual mode, with treatment time monitored at higher heater outputs.
- Radiant heater operated with automatic skin temperature control:  
Skin temperature controlled and monitored.

**Always maintain a minimum distance of 60 cm between the infant radiant heater and the bed. Otherwise, risk of burning the patient.**

**Do not use skin temperature control for patients in shock or with a high temperature.**

**Do not use mobile phones within 10 metres of the apparatus.**

Mobile phones can disrupt the function of electro-medical equipment and can thereby endanger the health or life of patients.\*

\* Dräger medical appliances comply with the interference immunity requirements of the specific standards for the products or EN 60601-1-2 (IEC 601-1-2). However, depending on the design of the mobile phone and situation of use, field strengths may occur in the immediate environment of a mobile phone that exceed the values of the standards quoted and therefore cause interference.

**Installation options:**

- Mounted on wheeled stand
- Wall-mounted

1-382-94

## Preparation

### Positioning the apparatus


#### Installation on wheeled stand

- 1 Unscrew the height-setting screw and set the distance between the radiant heater and the patient.  
**Always maintain the minimum clear height between the bottom edge of the infant radiant heater and the patient = 60 cm.**
- Retighten the height-setting screw. Check that the tube cannot slip down.
- Position the stand so that the radiant heat is directly above the patient.
- 2 Turn the swivel-adjusting screw and tilt the radiant heater so that the radiant heat is directed evenly at the patient.

#### Prevention of tipping

Do not tilt the apparatus by more than 10°.

- Do not place any objects on the trolley bracket.

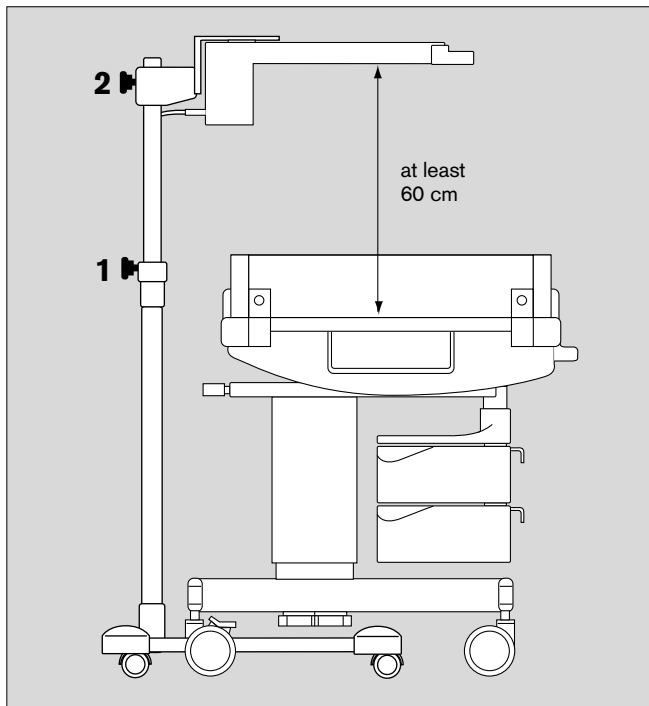
#### Wall-mounted installation

- Always keep a minimum clear distance of 60 cm to the patient.

### Checking readiness for operation

#### Before using for the first time

- Check that the mains power supply conforms to the specifications on the rating plate next to the power cord.
- Plug in the power plug.

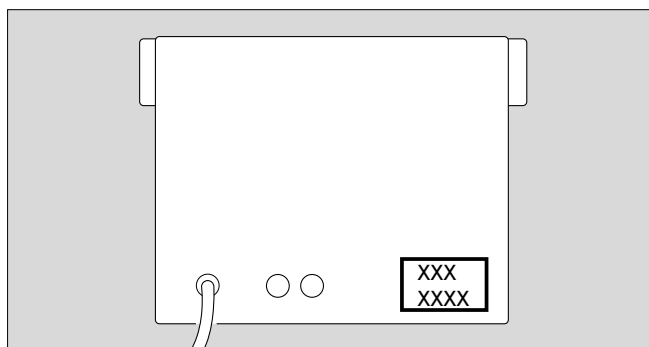
### Before each use


- 1 Press the power switch until it clicks into position = ON.

The microprocessor then runs a self-test to check that the apparatus is functioning correctly.  
All displays come on for about 1 second.  
All numeric displays show **88.8**.  
A continuous tone is sounded.  
The apparatus then switches over to »Man« mode (manual adjustment of heater output).

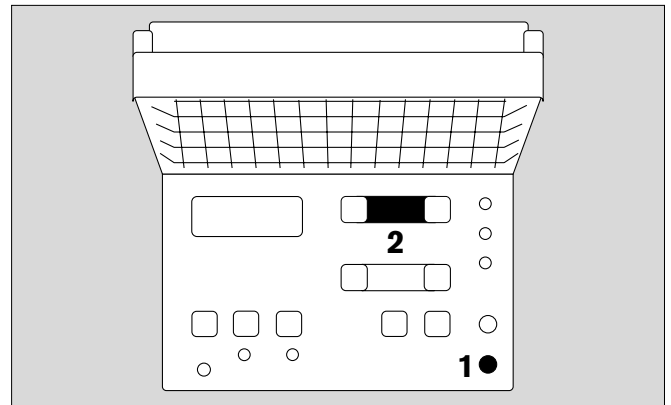
- 2 The message **Err** on the display indicates a fault. See page 50.

- 3 5 LEDs are lit in the heater output display = heater level 5.
- 4 Press key  or  and release immediately. Heater level 5 is accepted.  
If the keys are held down, the heater level is changed.
- 5 Press key : a reference temperature is simulated to check the temperature measurement.
- 6 The message **36 °C (±0.1 °C)** should appear in the temperature display. If not, see page 50.

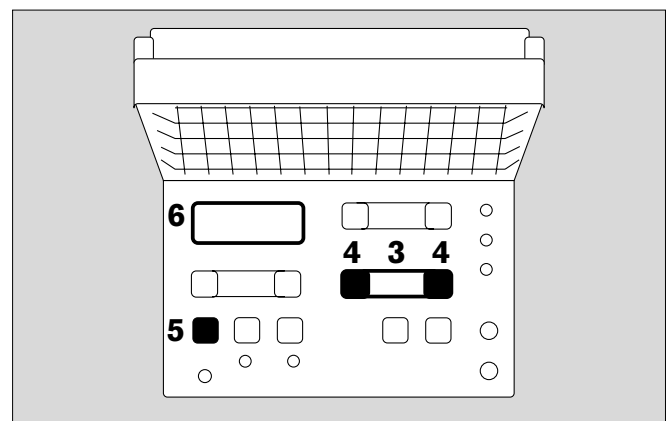
### Checking the power failure alarm and NiCd rechargeable battery

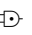- Disconnect the power plug.
- 7 The red LED marked  comes on. An intermittent tone is sounded: the volume of this signal must remain constant for at least 30 seconds.

If the volume fades earlier:

- Leave the apparatus switched on and connected to the mains until the battery is charged.
- Repeat the test.

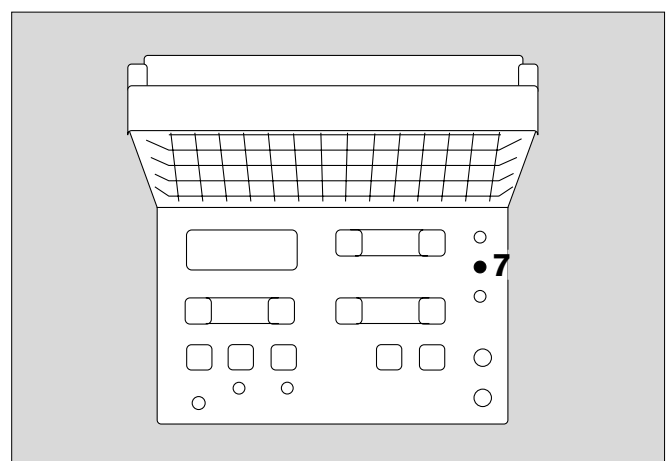
In the event of power failure, or if the apparatus is accidentally switched off, the settings remain stored in memory for 5 minutes.

## Preparation

### Checking readiness for operation

#### Checking the lighting

- 1 Press the light switch :  
The resting surface of the bed will be illuminated.
- 1 Press switch  once again:  
The light will be switched off.
- 2 Press the power switch in as far as it will go = OFF.
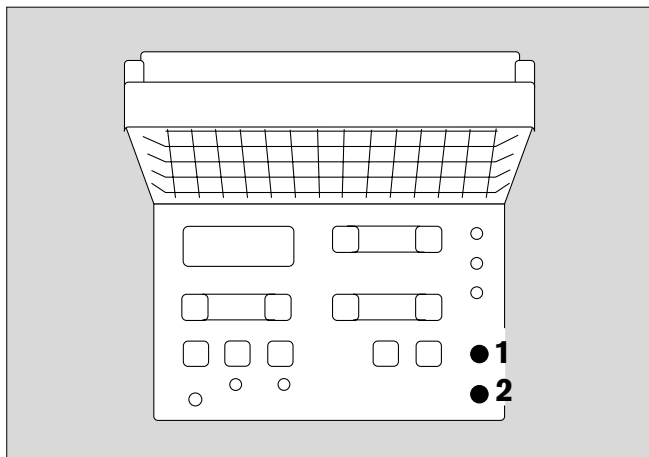
#### Checking the skin temperature sensor

Only for skin temperature control.

With the apparatus switched on,


- 3 Press the  key: the apparatus switches over to »Man« (manual) mode.
  - 4 Plug in the skin temperature sensor.
  - 5 The measured temperature is displayed,
- or

3 dashes at bottom of display = temperature less than 33 °C

3 dashes at top of display = temperature above 38 °C

If 3 dashes are displayed in the middle of the display:

- See section on "Troubleshooting", page 50.
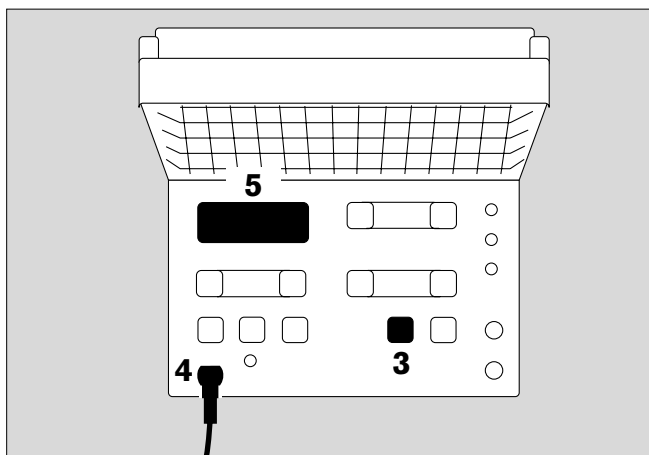
#### Checking that the care devices are ready for operation

- Check the Babytherm 8000 WB and Photo-Therapy 4000 as indicated in their specific Instructions for Use.

## Operation

### Safety precautions

**Do not place any articles, e.g. material, on top of the infant radiant heater.**

Ventilation would be impeded, and the infant radiant heater and deposited article could be damaged.

**Do not touch the top surface of the radiant heater or the mesh guard. Danger of burns.**

**During operation, do not use any inflammable cleaning agents or medicaments. Fire hazard.**

**Additional external heat sources, such as a radiant photo-therapy unit or mattress heater increase the energy supply to the patient and must be taken into account.**

**Environmental conditions, e.g. draughts, can affect the temperature balance of the patient.**

**Avoid direct sunlight and draughts.**

**Draughts reduce the effect of the infant radiant heater at the surface of the bed.**

**Always keep the minimum distance of 60 cm between the infant radiant heater and the bed. Otherwise there is a risk of causing burns to the patient.**

**Wait for the heater to warm up.**

The apparatus is only ready for use after 25 minutes.

**Do not place medicines and infusion solutions within the radiation range.**

**Do not use skin temperature control for patients in shock.**

The skin temperature of a patient in shock is much lower than normal.

The patients would be excessively warmed if the skin temperature control system were used. For patients in shock, adjust the heater output in »Man« mode and measure the core temperature every 15 minutes.

**Do not use skin temperature control for patients with a high temperature.**

The skin temperature will be much higher than normal. The patients would be overcooled by the skin temperature control.

For the connection of auxiliary devices to the patient, use only devices conforming to EN 60601-1.

The measuring functions of the apparatus must be inspected regularly. See the section on "Maintenance Intervals", page 49.

## Setting the heater output manually

In »Man« mode

E.g. for short-term treatment or for infants in shock and for whom automatic skin temperature control must not be used.

The radiant heater delivers a preset heat output, regardless of the body temperature of the infant.

**Do not operate the apparatus without medical supervision.**

**Regularly check the core temperature of the infant.**

To remind the user to monitor the core temperature constantly at higher heater output levels, an intermittent tone (3 short bleeps) is sounded every 15 minutes at heater level 6 and above.

Acknowledge this audible alarm by pressing the  key (no audible alarm is emitted in the lower range up to heater level 5).

- 1 Press the power switch until it clicks into place = ON.  
Wait for completion of the self-test.  
The apparatus then switches over to »Man« mode.
- 2 The green LED in the  key lights up.
- 3 5 green LEDs flash in the heat level display = heat level 5.  
This setting can be used as a standby setting to warm the resting surface of the bed.
- 4 Briefly press key  or  : the green LED segments will remain constantly lit.

### Wait for the warm-up time

Wait 25 minutes before using the apparatus.

- Warm up in »Man« mode and with heat level 5.

### Changing the setting



- 5 Hold down key  or . The heater level can be adjusted in discrete increments, i.e. segment by segment.  
Each LED segment represents a heater level of approx. 10 % of the maximum heater output.
- 6 The yellow LED  lights up or flashes according to the heating phases.

### Regularly check the core temperature

At heater level 6 and above, an intermittent audible alarm is emitted every 15 minutes to remind the operator to check the patient's core temperature. The radiant heater continues to heat.

- 7 Acknowledge this signal with the  key.

If a skin temperature sensor is connected or applied in manual mode, it will generate a skin temperature display but will not automatically control the skin temperature.
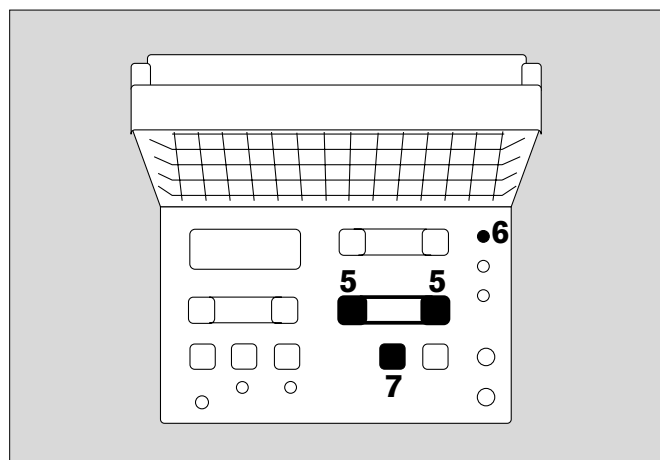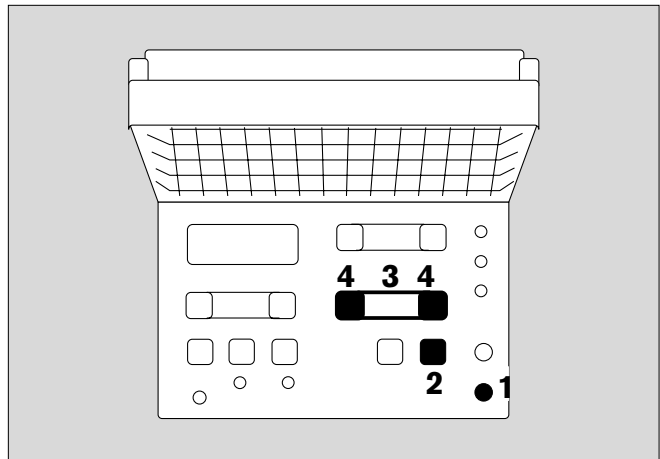
## Using skin temperature control

In this operating mode, the skin temperature of the infant is constantly adjusted towards a selected temperature setting. The sensor attached to the skin measures the skin temperature. If the measured skin temperature is lower than the set value, the radiant heater level is increased. If the skin temperature is higher than the set value, the radiant heater level is reduced.

Therefore:

**⚠ Do not use for infants in a state of shock. Their skin temperature will be much lower than normal. If automatic skin temperature control is used, it would cause overheating of the body (hyperthermia).**


Set the heater level manually, see page 38.

**⚠ Do not use for infants with a high temperature (fever). Their skin temperature will be much higher than normal. If automatic skin temperature control is used, it could induce hypothermia.**

Check the set value or set the heater level manually, see page 38.

When using the disc-shaped skin temperature sensor:


**1** Plug the sensor connector into the socket.

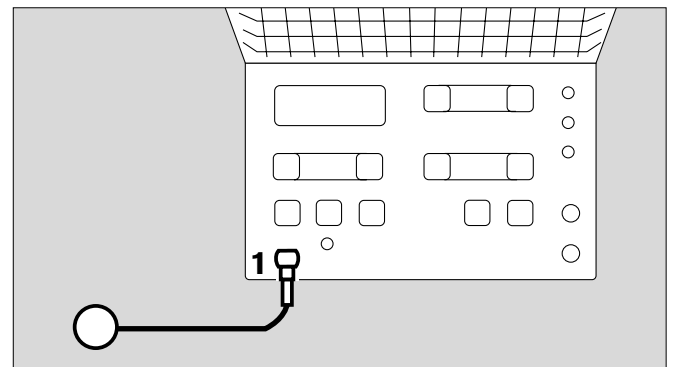When using the rod-shaped skin temperature sensor:

**2** Connect the adapter cable to the sensor cable and


**3** Plug the adapter connector into the socket.
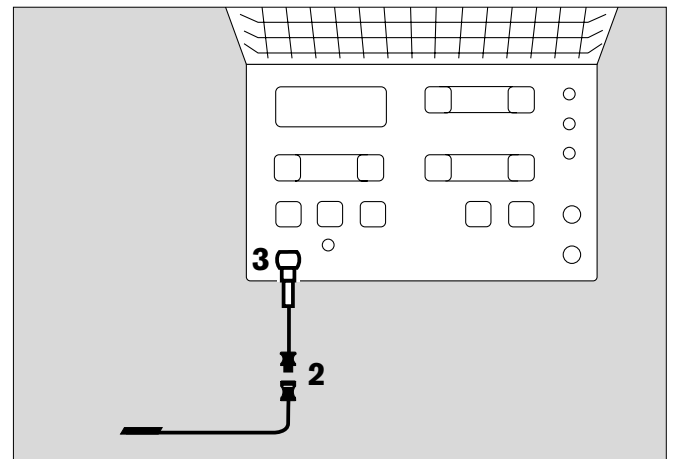
### Connecting the skin temperature sensor

If using the disc-shaped skin temperature sensor:

- Fit the double-sided adhesive ring to the temperature sensor:
- Detach the adhesive ring from its backing and stick it concentrically to the annular surface of the sensor.
- Detach the protective foil from the adhesive ring and stick the sensor to the appropriate part of the skin with the adhesive ring.

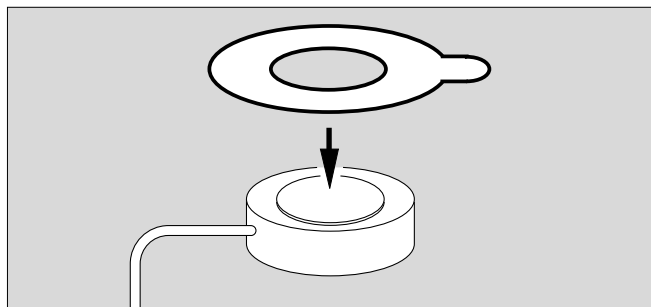If using the rod-shaped skin temperature sensor:

- Attach the tip of the sensor to the appropriate part of the skin with adhesive tape, e.g. 3M Microfoam 25 mm/1528-1.

If the infant is in the dorsal position:

- Attach the sensor to the abdomen in the region of the liver.


If the infant is in the prone position:

- Attach the sensor to the back, preferably in the region of the kidneys.

**Never place the sensor under the baby: otherwise, the core temperature rather than skin temperature would be measured and adjusted.**


- Secure the sensor cable with adhesive tape.
- **Regularly check that the skin temperature sensor is in the correct position.**  
If the skin temperature has fallen from the infant, it will measure air temperature, leading to a risk of excessively warming the infant.
- Never use the skin temperature sensor to measure rectal temperature.
- **During skin temperature control, the sensor must never be used rectally, since otherwise it would be measuring core temperature instead of skin temperature as feedback control parameter.**

- 1 Press the power switch until it clicks into position = ON.  
After the self-test, the apparatus switches to »Man« (manual) operating mode.

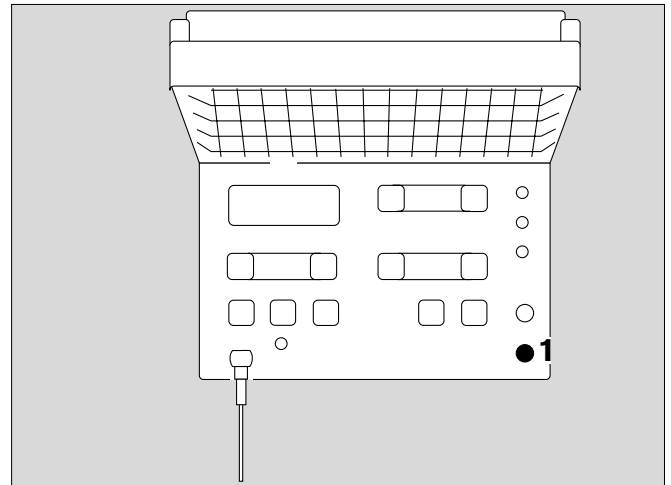If the temperature is outside the measuring range from 33 °C to 38 °C:


- 3 dashes at bottom of display = Temperature below 33 °C.
- 3 dashes at top of display = Temperature above 38 °C.
- See section on "Troubleshooting", page 50.

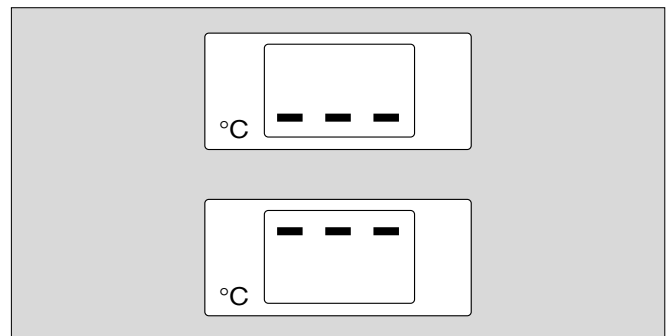**Allow at least 5 minutes for the skin temperature sensor to adjust to the temperature of the infant.**

If the displayed value remains constant:

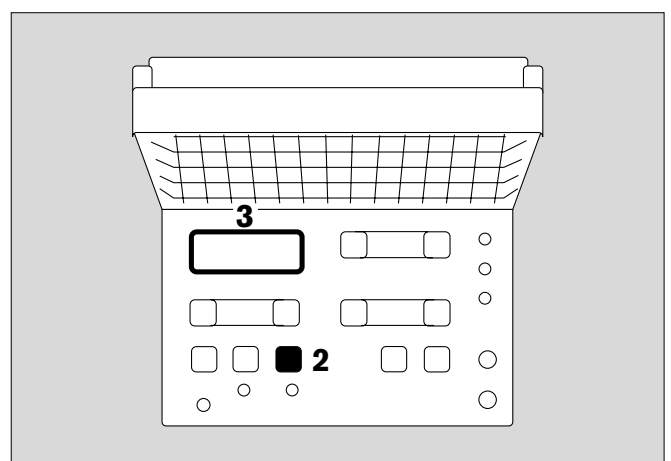- 2 Press the  key: skin temperature control is activated, and the green LED is lit.
- 3 The skin temperature is displayed in the temperature display.


## Operation


### Using skin temperature control

The system initially proposes the default temperature setting 36.3 °C.

- 1 The value **36.3 °C** flashes in the display.
- 2 Briefly press and release  or  to confirm the value.
- 1 The displayed value remains constant.


The temperature setting can be adjusted in 0.1 °C steps:

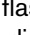- 2 Hold down key  or  until the desired temperature setting appears in the display.
- 3 The yellow LED  lights up or flashes according to the adjustment (heater function display).

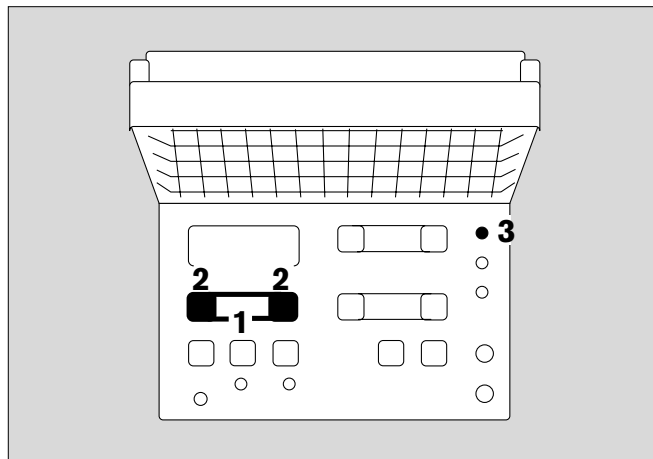To check heater function during operation:

- 4 Press the  key.
- 5 The simulated temperature value must be **36 °C** (±0.1 °C) (function check).

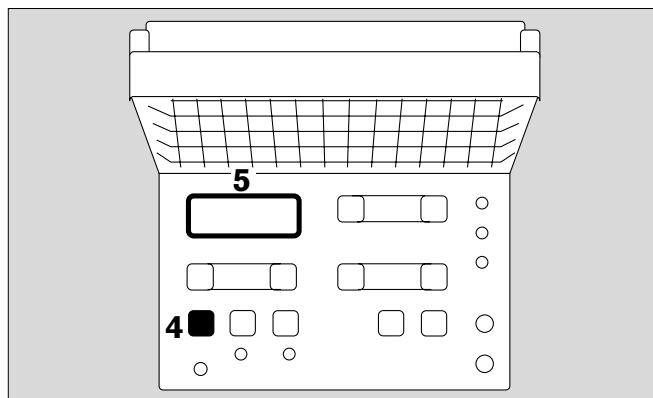**If the sensor probe is disconnected or the sensor is defective:**


- 6 Three flashing dashes in the middle of the display
- 7 the red alarm LED marked **Sensor** flashes.

A continuous audible alarm tone is emitted.

- Immediately connect the sensor plug or change the skin temperature sensor.

If the error cannot be remedied immediately:

- Remove the skin temperature sensor from the patient.
- Switch to manual mode (»Man«).
- See "Troubleshooting", page 50.
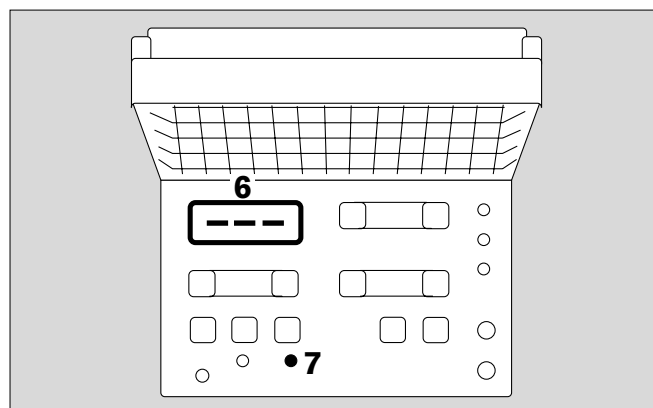
### Allow time for the system to reach steady state

Deviations between the set and measured skin temperature are normally **corrected within 5 to 15 minutes**.

An infant's skin temperature changes frequently, for instance as a result of food intake or handling. Deviations of a few tenths of a degree are normal.

Therefore:


Only change the set value for skin temperature if the core temperature needs to be corrected.

For deviations of more than  $\pm 0.5^{\circ}\text{C}$  between the set skin temperature and the measured skin temperature:

- 1 The red LED marked  $\pm 0.5^{\circ}\text{C}$  flashes, and an intermittent tone is sounded.

The intermittent tone can be muted for approx. 10 minutes:

- 2 Press the  key. The yellow LED in the key will come on.
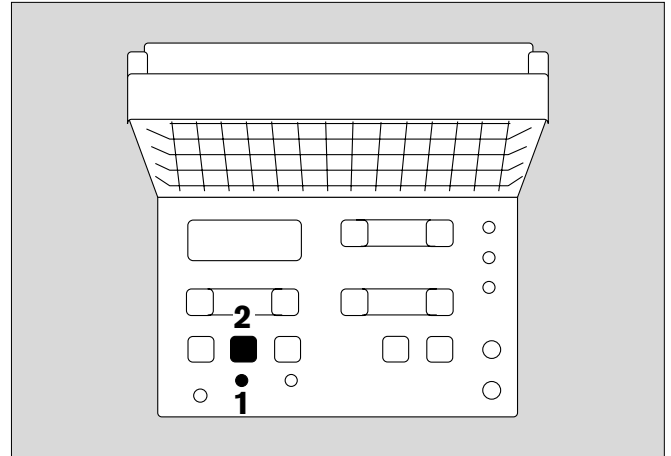
### Using the timer

for up to 100 minutes.

If the timer is started e.g. at the time of a birth, three short tones are sounded after 1 minute, after 5 minutes and after 10 minutes to prompt an entry in the APGAR report.

To start the timer:


- 3 Press key .

To stop:

- 4 Press key .

To reset the display:

- 4 Press key  a second time.

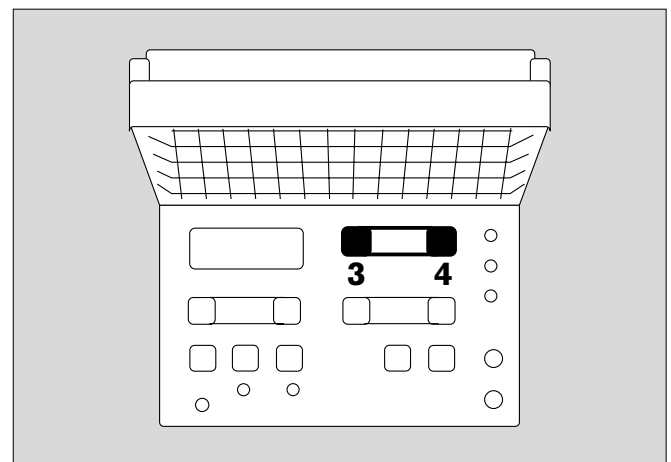
## Adjusting the halogen lamp beam

for precisely aimed lighting

- 1 Press the tabs towards the front panel to adjust the light cone towards the rear.

Pull the tabs

to move the light cone forwards.
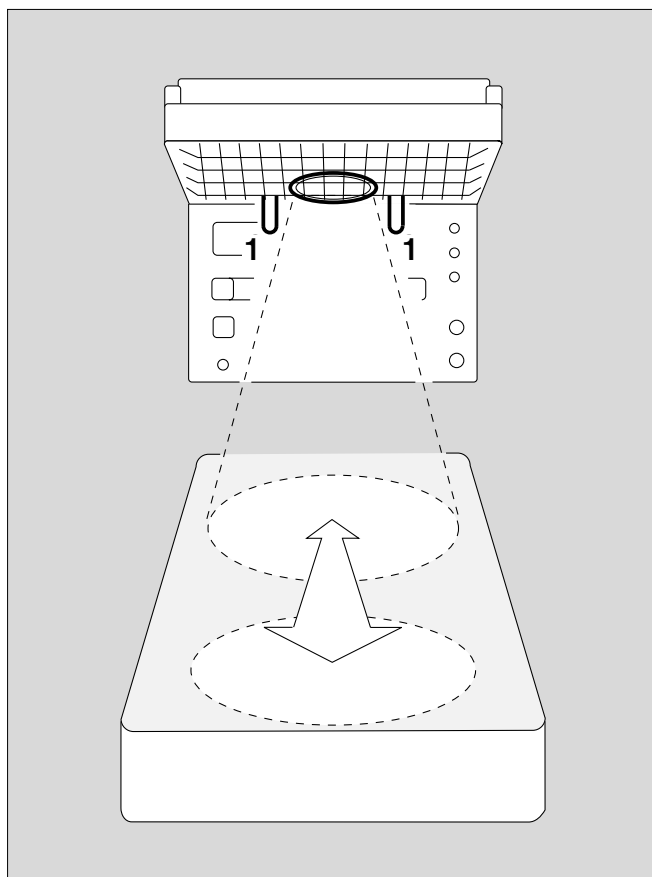
## Shutdown


First:

- Fully press in the power switch = OFF.

Then:

- Disconnect the power plug.  
These two steps must not be performed in reverse order, because otherwise the power failure alarm is activated.

Wait approx. 30 minutes for the apparatus to cool down before servicing/cleaning.

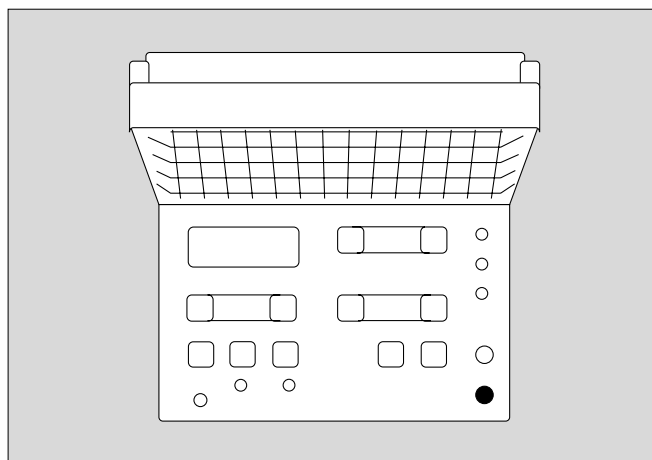
## Setting instructions


### Operating the RH 600 infant radiant heater

Constantly monitor the core temperature and adjust the temperature setting to the requirements of the patient.

- Do not fit a cover panel, otherwise the radiant heater will be ineffective.

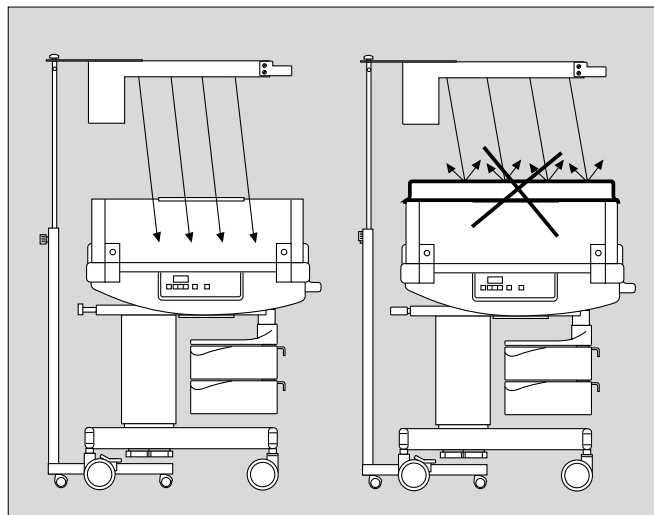
- Do not cover or clothe the infant.

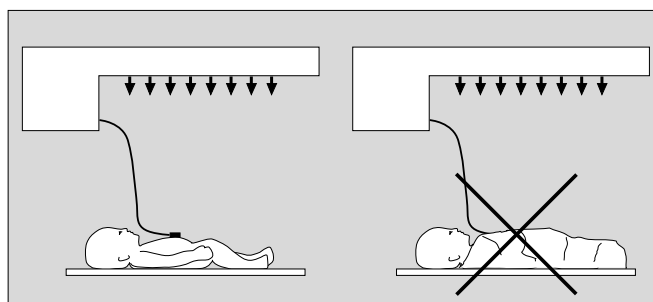
### Operation of the RH 600 infant radiant heater and mattress heater with gel mattress (e.g. Babytherm 8000 WB)

- Set the mattress temperature to the appropriate **core temperature** for the patient, e.g. between 38.0 °C and 38.3 °C for a premature neonate, or 37.0 °C for a neonate born after the full term.  
Set the skin temperature setting to the appropriate **skin temperature** for the patient, e.g. 37.0 °C for a premature neonate, or 36.5 °C for a neonate born full term.
- Increasing the temperature of the patient:  
Set the mattress temperature to the desired core temperature. If necessary, increase the skin temperature setting accordingly.


For infants with severe hypothermia: adjust the mattress temperature and skin temperature settings so that the patient is warmed by about 1 °C per hour.

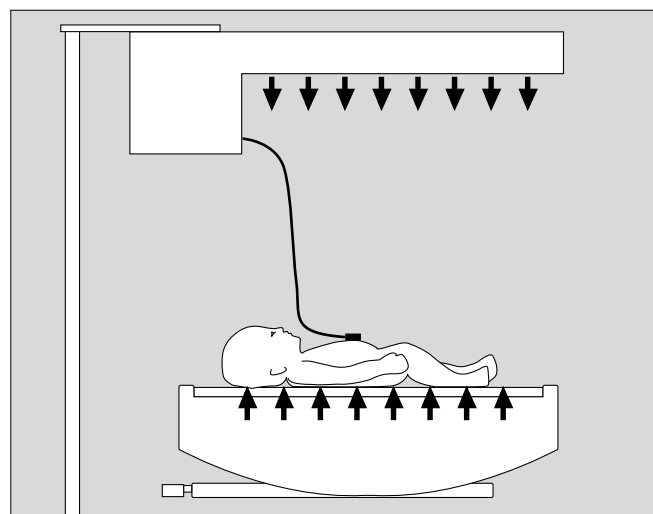- Stabilising the temperature of the patient:  
Set the mattress temperature to the present core temperature, and set the skin temperature control setting to the present skin temperature.

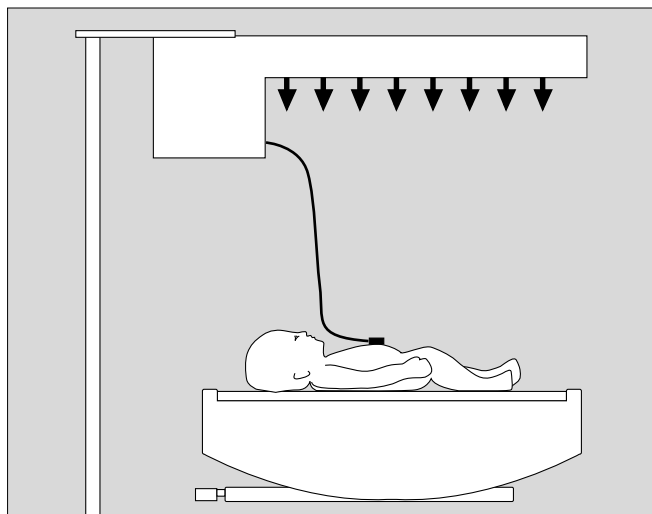- Reducing the temperature of the patient:  
Set the mattress temperature to the desired core temperature. If necessary, slightly reduce the skin temperature setting accordingly.
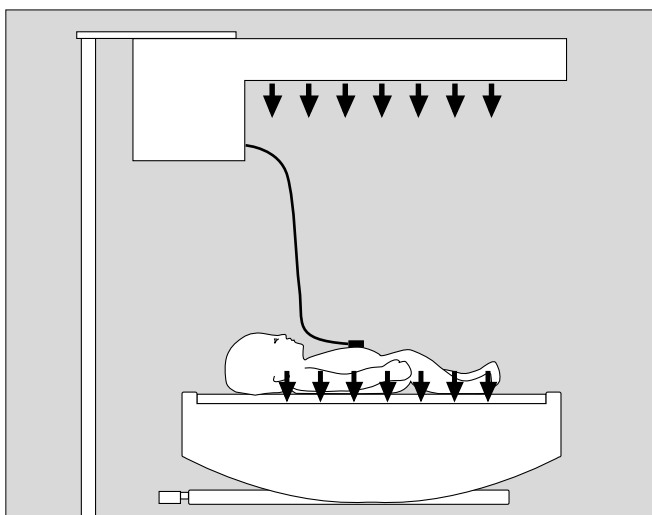
## Phototherapy

**To reduce risk of hyperthermie you may use phototherapy in combination with the RH 600 infant radiant heater.**

Otherwise, risk of inducing hypothermia in the undressed infant.

Example: Wall-mounted RH 600 with Babytherm 8000 WB and Photo-Therapy 4000


- Push the stand of the phototherapy unit under the frame of the infant bed, e.g. Babytherm 8000 WB.
- Tilt the lamp housing towards the bed.
- Adjust the height of the phototherapy unit: **minimum clearance** between the lower edge of the lamp body and the patient: **0.3 metres**.
- Position the lamp housing so that the radiant heat reaches the patient unimpeded.
- Follow the specific Instructions for Use of the phototherapy units.

If operating the radiant heater in »Man« mode:

- Reduce the heater output level by 3 steps compared to operation without phototherapy.

Using the phototherapy units will result in increasing the temperature of the infant.

**Monitor the core temperature of the infant.**
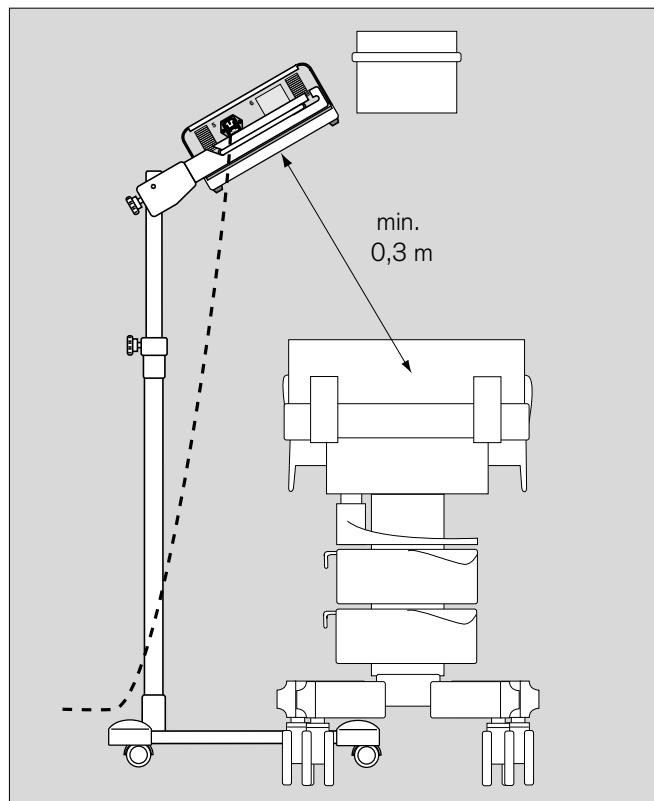
## Care

### Shutdown

- Switch off the device(s). Disconnect the power plug.

### Cleaning/disinfecting

- Allow the radiant heater to cool down for about 30 minutes before cleaning/disinfecting.
- Do not allow any liquids to enter the apparatus.

#### Wheeled stand;

#### Infant radiant heater, exterior:

- Remove any visible impurities with a disposable cloth soaked in washing-up liquid.
- Wipe-disinfect all surfaces.
- Wait the specified time for the detergent to act and then wipe off with a damp cloth and dry.

Only use the recommended detergents and disinfectants. Unauthorised products, e.g. alcohol, may cause stress cracking of the acrylic glass and "macrolon" material.

Use only preparations classified as "surface disinfectants" for disinfecting. For material compatibility, we recommend preparations based on

- aldehydes
- quaternary ammonium compounds.

Due to possible damage to materials, the following preparations are unsuitable:

- halogen-releasing compounds
- strong organic acids
- oxygen-releasing compounds.

For users in the Federal Republic of Germany, we recommend the use of disinfectants listed in the current DGHM list (DGHM = German Society for Hygiene and Microbiology).

The DGHM list (published by mhp Verlag GmbH, Wiesbaden) also specifies the active basis of each disinfectant.

For countries where the DGHM list is unavailable, we recommend products based on the above active compounds.

For example, the following wipe-disinfectants may be used:

|                 |                                    |
|-----------------|------------------------------------|
| Bacillotox      | Bode, Hamburg                      |
| Incidur         | Henkel Hygiene GmbH,<br>Düsseldorf |
| Dismozon pur    | Bode Chemie GmbH & Co.             |
| Sekusept-Pulver | Henkel Hygiene GmbH,<br>Düsseldorf |

Follow the Instructions for Use provided by the detergent manufacturer.

### Before next use

- After wipe-disinfection, operate the fully assembled unit for a few hours without the patient, in order to eliminate any disinfectant residues.  
Set the radiant heater to heat level 5 in »Man« mode.
- Check that the unit is ready for operation:  
see page 34.

## Maintenance intervals

**Always clean and disinfect the unit and accessories before any maintenance – even when returning the unit to supplier for repair.**

**⚠ Always disconnect power supply before any maintenance.**

|                                       |   |
|---------------------------------------|---|
| Temperature measuring system          | The complete measuring chain must be checked as part of the safety check by factory-trained technicians |
| Ni-Cd battery for power failure alarm | Must be replaced after 2 years by factory-trained technicians   |
| Time Keeper RAM                       | After 4 years by factory-trained technicians  |
| Metrological inspections              | Once a year by factory-trained technicians  |
| Inspection and maintenance            | Once a year by factory-trained technicians  |

## Disposal of used batteries

- Do not throw into fire.
- Do not force open.

Danger of chemical burns. Batteries may explode.

Batteries must be treated as special waste:


- Discard in conformity with local waste disposal regulations.

Please obtain information from your local environmental and planning authorities and waste disposal companies.

## Troubleshooting

| Fault indicator  | Cause  | Remedy  |
|--|--|---|
| Red LED marked $\pm 0,5\text{ }^{\circ}\text{C}$ flashes.<br>Intermittent tone is sounded.   | Actual value of skin temperature deviates by more than $\pm 0.5\text{ }^{\circ}\text{C}$ from the set value. | If the measured value is less than the set value:<br>check that the skin temperature sensor is correctly fixed.<br><br>If the measured value is greater than the set value:<br><b>measure the core temperature.</b> |
| Red LED marked <b>Sensor</b> comes on.<br>Continuous tone.<br>In the skin temperature display:<br>3 dashes in the middle.  | Connector of the skin temperature sensor not plugged in.<br><br>Sensor defective.                            | Check the plug-in connection.<br><br>Remove the sensor from the patient; switch over to »Man« mode; replace sensor, and then switch back to »Auto« mode.  |
| Red LED marked  comes on.<br>Continuous tone.   | No electric voltage.   | Check that the power plug is connected.<br>Check that both mains fuses are working correctly.<br>Check that the mains voltage is available.<br>Call Dräger Service.   |
| Red LED marked <b>Inop</b> comes on.<br>In the timer, the message <b>Err</b> is displayed alternately with an error code.<br>Constant tone starts.   | Malfunction  | Note the error code.<br>Call Dräger Service.  |
| Lighting does not work.  | Halogen lamp defective.  | Have the halogen lamp replaced by trained personnel.  |
| Insufficient radiant power.  | One or more infrared radiators defective.  | Have the radiant heater replaced by trained personnel.  |
| 36 °C reference temperature is outside the tolerance of $\pm 0.1\text{ }^{\circ}\text{C}$ when key  is pressed. | Inaccurate operation of the temperature measuring device.  | Call Dräger Service.  |
| In the skin temperature display:<br>3 dashes at bottom<br>3 dashes at top  | Skin temperature $< 33\text{ }^{\circ}\text{C}$<br>Skin temperature $> 38\text{ }^{\circ}\text{C}$ .         | Check that the skin temperature sensor is correctly fixed.  |

## What's what

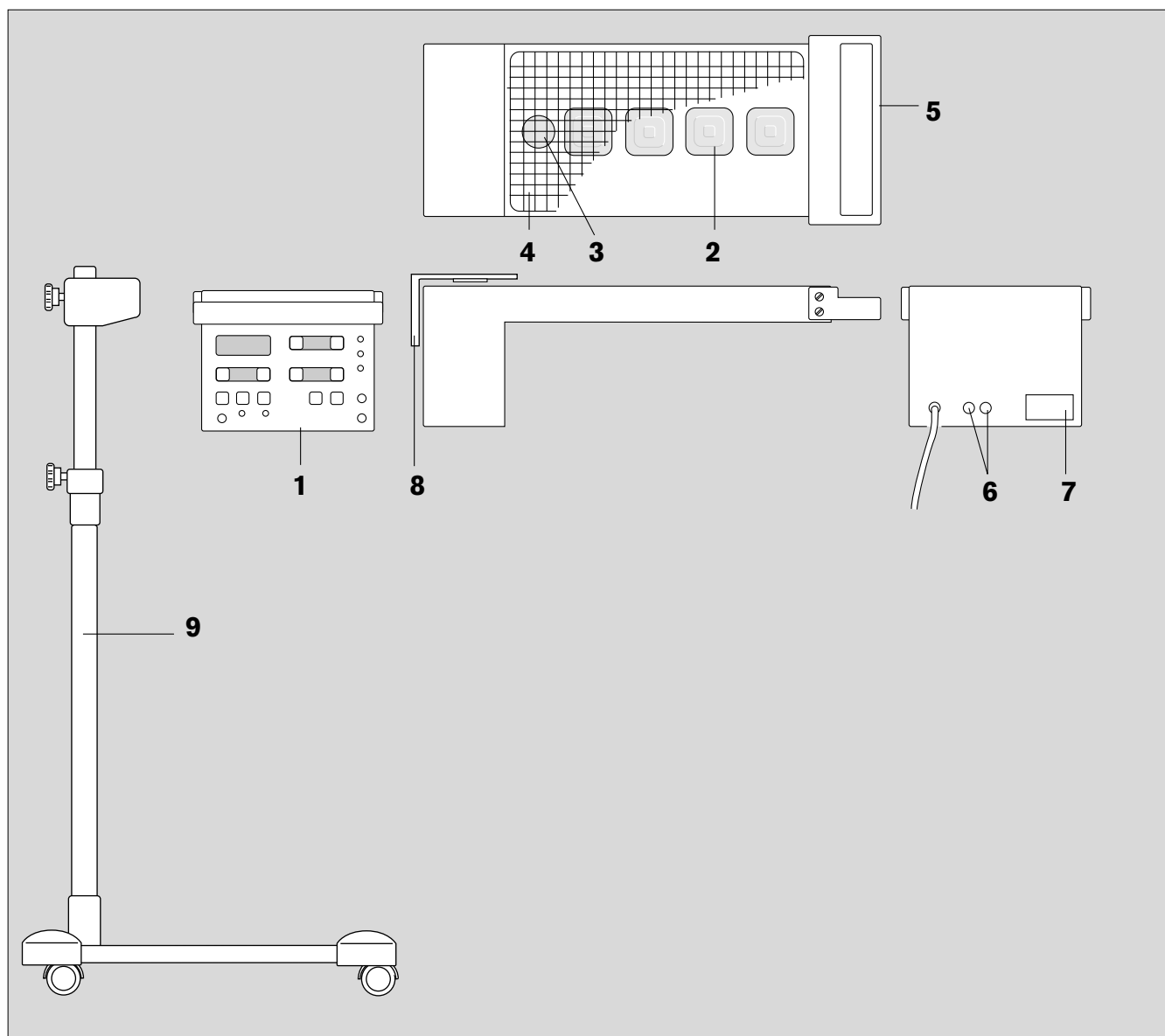
- 1 Control panel of the radiant heater
- 2 Infrared ceramic radiating elements (4x)
- 3 Halogen lamp
- 4 Protective mesh
- 5 Handle
- 6 System fuse (2x)
- 7 Rating plate
- 8 Wall bracket

- 9 Wheeled stand with fixing screw for height adjustment and setscrew for setting the tilt angle of the radiant heater housing.

⚠ Follow instructions.

## Control panel

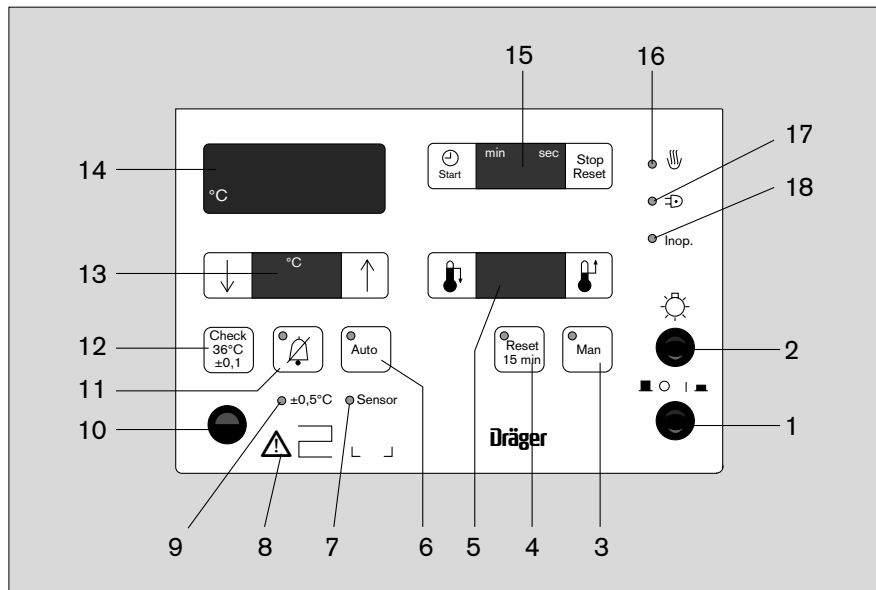- 1 Power switch
- 2 Lamp switch
- 3 Push-button for »Man« (manual) mode;  
when activated, the green LED is lit.
- 4 Push-button for acknowledging the intermittent  
alarm tone that is sounded every 15 minutes at  
heater level 5 and above to remind the user to  
measure the core temperature.  
The yellow LED goes out after the alarm is  
acknowledged.
- 5 Green bar-segment display for the heater output  
level; one segment corresponds to a heater level  
of 10% max. heater power;  
To the left of the bar display:  
Key for reducing the heater level.  
To the right of the bar display:  
Key for increasing the heater level.
- 6 Push-button for »Auto« (automatic) mode;  
the green LED comes on when the button is  
pressed = skin temperature control activated.
- 7 Red LED "Sensor":  
comes on during skin temperature control if the  
skin temperature sensor is removed or defective.
- 8 Follow the Instructions for Use of the skin  
temperature control.
- 9 Red LED  $\pm 0,5^{\circ}\text{C}$ ;  
flashes when the measured skin temperature  
deviates by more than  $\pm 0,5^{\circ}\text{C}$  from the set value.
- 10 Connection for skin temperature sensor
- 11 Push-button for muting the intermittent alarm tone for  
10 minutes;  
the yellow LED comes on while the intermittent tone  
is muted.
- 12 Push-button for simulating the temperature  
 $36 \pm 0,1^{\circ}\text{C}$  (function check).
- 13 Display for the set value of the skin temperature.  
To the left of setpoint display:  
Key for reducing the skin temperature setting.  
To the right of setpoint display:  
Key for increasing the skin temperature setting.
- 14 Display for the actual measured skin temperature.
- 15 Timer display.  
To the left of the display:  
Start key.  
To the right of the display:  
Stop key (1x)  
Reset key (2x)
- 16 Yellow LED ;  
lights up/flashes according to heater phases.
- 17 Red LED ;  
lights up in the event of a power failure.
- 18 Red LED Inop;  
lights up in the event of a malfunction.

## Technical Data

### Ambient conditions

During operation:

|                      |                            |
|----------------------|----------------------------|
| Temperature          | 15 °C to 30 °C             |
| Atmospheric pressure | 900 hPa to 1100 hPa        |
| Relative humidity    | 0 to 90 %, no condensation |
| Air velocity         | Max. 0.3 m/s               |

During storage/transport:

|                      |                            |
|----------------------|----------------------------|
| Temperature          | –20 °C to +60 °C           |
| Atmospheric pressure | 900 hPa to 1100 hPa        |
| Relative humidity    | 0 to 90 %, no condensation |

### Radiant heater

Radiant power

| Distance between radiant heater and resting surface of bed, in cm | Radiated power at heater level 5, in mW/cm <sup>2</sup> | Radiated power at heater level 10, in mW/cm <sup>2</sup> |
|---|---|--|
| 80  | 7   | 14   |
| 70  | 9   | 18   |
| 60  | 10  | 25   |

Light

Halogen, 12 V / 35 W

Timer

0 to 99 minutes 59 seconds;  
Start and stop times programmed in seconds;  
Reset function;  
Audible signal at 1 minute, 5 minutes and 10 minutes

Clearance between top edge of radiant heater housing and ceiling

>50 cm

### Skin temperature measurement

|   |  |
|---|--|
| Sensor                                    | Use only genuine Dräger earthed sensors:<br>see Order List on page 57. |
| Measuring range                           | 33 °C to 38 °C   |
| Precision                                 | ±0.1 °C  |
| Precision of the complete measuring chain | ±0.2 °C  |
| Setting range                             | 35 °C to 37 °C   |

### Operating characteristics

|                         |  |
|-------------------------|--|
| Mains voltage           | 100/120/127/230/240 V (please specify on order)<br>50/60 Hz                                |
| Power input             | 680 W  |
| System fuses            | M 6,3 A L 250 V; IEC 127-2/III (2 fuses)   |
| Protection class        | I  |
| Skin temperature sensor |  Type B |
| Standards               | EN 60601-1<br>EN 60601-2-21  |


### Electromagnetic compatibility EMC

Tested to EN 60601-1-2


|                                  |      |
|----------------------------------|------|
| Weight<br>(without fixing clamp) | 7 kg |
|----------------------------------|------|

## Dimensions

### Radiant heater
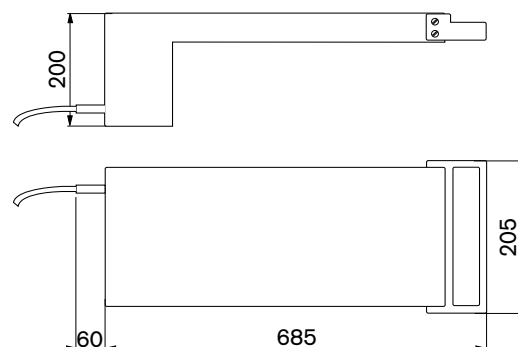
### Radiant heater on wheeled stand
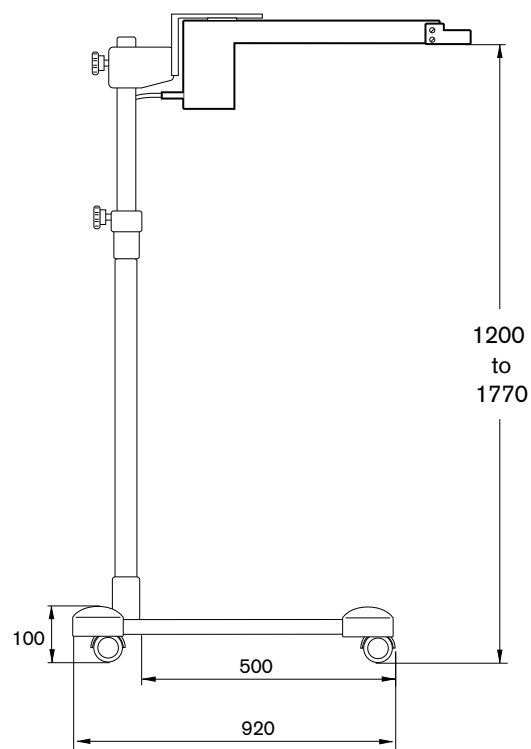
## Construction and description

The infant radiant heater consists of the heater unit and the control unit.

The heater unit comprises 4 infrared ceramic radiating elements and a low-voltage halogen lamp.

The halogen light provides good illumination of the working area and, compared to incandescent lights, good colour reproduction, so that diagnostically significant colour shades can be more easily recognised.

The control unit is located on the side of the housing. The part of the control unit facing the work area has a membrane keypad with the display and input keys.

In manual control mode, the heater level is set to fixed values.

With skin temperature control (automatic mode), the skin temperature is continuously adjusted towards a predefined set value. In both manual and automatic modes, correct monitoring is ensured by control routines and alarms.

With the built-in timer, started at the time of birth, the times relevant for the APGAR test (1 minute, 5 minutes and 10 minutes) are indicated by an audible signal.

## Order List

| Name and description   | Part No.         |
|--|------------------|
| <b>RH 600 Infant Radiant Heater</b><br>on wheeled stand<br>Mains voltage to order<br>100, 120, 127, 230 or 240 V; 50/60 Hz<br>with manual heater output control,<br>Skin temperature control,<br>Halogen light and timer | <b>2M 20 539</b> |
| <b>RH 600 Infant Radiant Heater</b><br>for wall-mounting   | <b>2M 20 116</b> |
| Set of fixings for wall mounting   | 2M 20 540        |
| <b>Consumables</b>   |                  |
| Skin temperature sensor,<br>incl. 250 adhesive rings   | 82 00 757        |
| Adhesive rings, set of 250   | 68 50 021        |
| Skin temperature sensor, yellow, set of 10   | 2M 21 916        |
| Adhesive pads, set of 20   | 2M 21 735        |
| Replacement adhesive pads, set of 100  | 2M 21 734        |
| Accessories required for rod-type<br>skin temperature sensor 2M 21 916:<br>Adapter cable   | 2M 20 736        |
| Halogen light, 12 V, 35 W  | 18 33 537        |
| <b>Phototherapy</b>  |                  |
| Photo-Therapy 4000, 230 V  | 2M 21 000        |
| Photo-Therapy 4000, 120 V  | 2M 21 700        |
| Stand for phototherapy unit  | 2M 18 780        |
| <b>Warming bed</b>   |                  |
| Babytherm 8000 WB 230 V  | FR 00 050        |
| Babytherm 8000 WB 120 V  | FR 00 088        |
| Technical documentation available<br>on request  |                  |

## Parts List

As alternatives to the part numbers listed in the "Order List", the following parts or devices that are no longer included in the standard Dräger equipment and accessories may be used:

| Name and Description                     | Order No. |
|--|-----------|
| PT 8000                                  | 2M 20 520 |
| PT 800                                   | 2M 18 800 |
| Phototherapy stand                       | 2M 18 780 |
| Skin temperature sensor, Export          | 82 00 757 |
| Skin temperature sensor, calibrated      | 82 00 756 |
| Set, skin temperature sensor             | 2M 21 726 |
| Set, skin temperature sensor, calibrated | 2M 22 010 |

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